

MYTREX

Distributor: SOTSU MEDICAL Co., Ltd.
3-1-10-5F, Hakataekimammi, Hakata-ku, Fukuoka City, Fukuoka
812-0016, JAPAN
https://mytrex.jp

©2025 SOTSU MEDICAL CO., LTD.
20250908B

MYTREX RAKUNO FIT

Instruction Manual

For Storage
Includes Warranty Card



Thank you very much for purchasing this product. Please read this instruction manual carefully and use the product correctly and safely. This instruction manual includes a warranty card. Please keep it properly even after reading it.

Safety Precautions

Please follow the prescribed procedures and use the product without overexertion. If any abnormalities occur, immediately stop using the product and contact our customer support. Continuing to use the product in a malfunctioning state may pose a risk of adverse health effects and must be strictly avoided. This section explains important safety precautions that must be strictly followed to prevent injury to the user or others, as well as damage to property.

About Warning Labels

The following symbols are used to indicate the severity of hazards or damage that may result from improper use.

- WARNING** Indicates "a risk of death or serious injury" if the product is used incorrectly.
- CAUTION** Indicates "a risk of minor injury or property damage" if the product is used incorrectly.
- Indicates actions that are strictly "prohibited".
- Indicates "mandatory" actions that must be strictly followed.

WARNING

Do not use the product while wearing accessories or other hard objects in the area of use.
This may cause accidents or injuries.

Do not use the product while driving a car or bicycle, or while performing other tasks. Do not stand or walk while using it.
This may cause accidents or injuries.

Do not use the product under a blanket. Do not use it while smoking or near heating appliances or on heated carpets.
This may cause accidents, electric shock, or low-temperature burns.

Do not allow pins, foreign objects, dust, or moisture to adhere to the connector or gaps of the main unit or accessories.
This may cause electric shock, short circuits, or fire.

MYTREX Customer Support

For questions or concerns about the product, please feel free to contact us.

Business hours: 9:00 - 17:30
(excluding Saturdays, Sundays, national holidays, summer holidays, and year-end/New Year holidays)

[Contact by phone]
092-260-7366

[Contact via the web]
https://mytrex.jp/support/



WARNING

Do not use in the following cases.
This may cause accidents or injuries.

- Individuals who have been prohibited from receiving massage by a doctor.
- Examples: thrombosis, embolism, severe aneurysms, acute varicose veins, various forms of dermatitis, skin infections (including subcutaneous inflammation), etc.
- Individuals with severe circulatory disorders in the legs.

Do not allow by individuals who are unable to express their intentions or operate the product.
This may cause accidents or injuries.

Do not allow small children to use or play with this product. Do not use on pets.
This may cause accidents, injuries, or health issues.

Do not use on any body part other than the calf. Do not use on the neck.
This may cause accidents or injuries.

WARNING

Do not use a deteriorated lithium-ion battery.
This may cause electric shock, short circuit, or fire.

Do not apply excessive force to the USB cable, such as damaging, modifying, or trying to fix it up. Also, do not place heavy objects on it, bend it, or pull it.
This may cause electric shock, short circuit, or fire.

When applying the USB cable, hold the connector and do not pull the cable. Also, do not unplug with wet hands.
This may cause electric shock, short circuit, or fire.

WARNING

Since this product contains a lithium-ion battery, the following actions are strictly prohibited.
This may cause accidents or injuries.

Do not use in the following cases.
This may cause accidents or injuries.

- Applying heat or throwing into fire
- Dropping, crushing, or exposing to moisture
- Using, charging, or leaving the product in high or low temperature environment
- Using, charging, or leaving the product on a blanket or near heating appliances or on heated carpets
- Using, charging, or leaving the product in high-temperature areas, such as in a car

If you notice any of the following abnormalities, stop using the product immediately and contact customer support.
This may cause accidents or injuries.

- The battery is swollen, making strange noises, or emitting unusual smells.
- The product does not operate even after charging.
- The product's operating time becomes extremely short.
- The product becomes excessively hot during charging.

Warranty Card

This warranty card guarantees free repair or replacement in the event of a malfunction if the product is used according to the instruction manual. If a malfunction occurs during the warranty period from the date of purchase, please present the product, warranty card, and a certificate of purchase receipt, delivery note, etc.) to the retailers or our customer support. For repairs after the warranty period, the function can be maintained by repairing the product, we will repair it for the fee the customer's request.

1. Even during the warranty period, this warranty does not apply to accidents caused by the following reasons or similar situations:
- Malfunction or damage caused by improper usage or storage, contrary to the instruction manual.
- Malfunction or damage due to improper handling (misuse, dropping, etc.), modification, or repair.
- Malfunction or damage caused by natural disasters (rain-related such as fire or lightning, environmental pollution) and damage, gas damage, abnormal voltage, or other external factors.
- Replacement of consumable parts such as batteries, due to wear and tear.
- Malfunction or damage when used outside of a general household setting.
- Malfunction or damage when purchased from unauthorized resellers or distributors.

2. Even during the warranty period, the customer is generally responsible for shipping costs when sending the product to the retailer.
3. When returning the product, please send it back in the same packaging conditions as received, including parts and packaging.

4. This warranty card cannot be returned, so please keep it in a safe place.
5. Our company reserves the right to amend and make the final interpretation of the warranty terms.

If the product malfunctions, please contact our customer support.
This warranty card assures the repair or replacement under the specified period and conditions and does not limit your rights.
*Minimum storage period for performance-related spare parts for this product will be limited for a maximum of 3 years after the product is disposed. These parts are necessary to maintain the functionality of the product.

Warranty Card
Model Numbers: MY-60K-2500 (Black), MY-60K-25P (Pink), MY-60K-25B (Blue)
Warranty period is 1 year from the date of purchase. (for the main unit)

Customer's Name/Address/TEL

Dealer's Name/Address/TEL

Purchase Date: Year: _____ Month: _____ Day: _____

WARNING

Consult a doctor before use if you fall under any of the following categories.
Using the product may cause accidents, injuries, or health issues. Even if none of the following applies, consult a doctor if you have any concerns about using the product.

● Individuals using implanted medical electrical devices that are susceptible to electromagnetic interference, such as pacemakers.
● Individuals with malignant tumors.
● Individuals with heart disease, heart conditions, or thrombosis.
● Individuals with circulatory disorders, respiratory dilation, or abnormal blood pressure.
● Individuals who are pregnant, or have recently given birth.
● Individuals currently taking medication.
● Individuals with sensory impairment due to advanced peripheral circulatory disorders such as diabetes.
● Individuals with skin abnormalities such as wounds, eczema, or rashes.
● Individuals with impaired temperature sensation.
● Individuals requiring braces or those with disabilities in their limbs.

WARNING

Do not intentionally remove, tear, or cut the fabric part of the main unit.
This may cause accidents, electric shock, or injury.

Do not splash water on the main unit or accessories, or use them in humid places such as bathrooms or outdoors.
This may cause fire, electric shock, or injury.

Do not disassemble, repair, or modify the product.
This may cause fire, malfunction, or accidents.

Materials
ABS resin, Polyester, Nylon, PP, Silicone, etc.
Detachable cover: 100% Polyester

Compatible Size
Calf circumference: Up to 48cm

Operating Environment
Temperature: 5°C to 40°C
Relative Humidity: 25% RH

Rated Time
10 minutes

Manufacturer/Distributor
SOTSU MEDICAL Co., Ltd.

Country of Manufacture
China

*Measured at 25°C under no-load operation. Actual time may vary depending on the usage environment and conditions.

Product Specifications

Brand	MYTREX
Product Name	RAKUNO FIT
Model Number	MY-60K-2500 (Black), MY-60K-25P (Pink), MY-60K-25B (Blue)
Rated Power Supply	5V=1A (Only with Class II symbol (□))
Rated Power Consumption	9.6W (Heat: 5W)
Battery	Lithium-ion battery 3.7V=4500mAh
Size	Approx. 135mm x 62mm x 45mm
Weight	Approx. 650g
Operating Duration*	Maximum operating time: Approx. 5.5 hours Minimum operating time: Approx. 3.5 hours

WARNING

Consult a doctor before use if you fall under any of the following categories.
Using the product may cause accidents, injuries, or health issues. Even if none of the following applies, consult a doctor if you have any concerns about using the product.

● Individuals with body temperature of 38°C or higher (feverish period) (e.g. those experiencing strong acute inflammatory symptoms such as fatigue, chills, or blood pressure fluctuations), or who are physically weak.

● Individuals who are feeling unwell or have abnormal physical conditions.
● Individuals with acute conditions (pain-related) such as fractures, osteoporosis, spinal fractures, sprains, or muscle strains.
● Individuals with spinal abnormalities or scoliosis.
● Individuals currently undergoing treatment as a medical institution.
● Individuals using the product on areas previously treated, operated on, or injected with drugs.
● Individuals who are extremely physically weak.
● Individuals with cold sores, herpes, or sores who use the product.
● Individuals with nerve disorders, skin sensitivity disorders, allergic constitution, skin diseases, atopic dermatitis, or sensitive skin.

WARNING

Immediately stop using the product and consult a doctor if you feel any abnormality during use.
Continued use may be harmful to your health.

Stop using the product if the main unit or accessories emit a strange odor or abnormal heat.
This may cause damage or fire.

Check for any tears in the detachable cover before use. Do not use if there is any damage.
This may cause accidents or injuries.

Use only the included USB cable.
Using a cable not included with the product may cause a malfunction. Do not use the accessories with other products.

Problems **Possible cause and solutions**

The product is making an unusual noise or a loud sound.
- The operating sound is generated as part of the product's design and is not considered abnormal.

The product cannot be charged.
- Are the AC adapter (sold separately) and USB cable connected correctly?
- Check the connection method.
- Charging may not be possible in high or low temperature environments.
- Charge the product in the recommended temperature range.

Troubleshooting

Before contacting us, please check the following again. If the product still does not function properly, please contact the retailer where you purchased the product or our customer support. Do not disassemble, repair, or modify the product under any circumstances.

Problems	Possible cause and solutions
The product does not operate.	- Is the battery level low? → Charge the product. - Have you been using it continuously? → Turn off the power and let the temperature return to normal. Individuals who have been using the product for a long time may overheat and trigger the protection function. Wait for a while before using it. - Is the battery level low? → Charge the product. - If approximately 10 minutes have passed since pressing the power button, the auto-off timer may have activated. This is not a malfunction.
The product suddenly stops.	- Is the battery level low? → Charge the product. - If approximately 10 minutes have passed since pressing the power button, the auto-off timer may have activated. This is not a malfunction.

WARNING

If the product malfunctions, immediately stop using it and turn off the power. If the product malfunctions during charging, or if there is a thunder or a power outage, unplug the AC adapter from the outlet immediately.
Do not continue to use the product if it is in an abnormal condition, such as emitting smoke, strange noises, or odors. This may cause accidents or fire.

CAUTION

Do not use the product for any purpose other than its intended use.
This may cause injury or health issues.

Do not use the product in the following situations.
This may cause accidents, injuries, or health issues.
● When extremely fatigued ● After consuming alcohol
● When tired before or after eating ● When the skin is abnormal

Do not charge, use, or store the product within the reach of infants or children.
This may cause accidents or injuries.

Do not insert hands or fingers into the moving parts during operation.
This may cause injury or health issues.

This product is designed for one foot only. Do not insert both feet.
This may cause accidents or injuries.

The main unit contains a built-in lithium-ion battery. When disposing of the product, please follow your local municipality's waste separation rules. For detailed regulations, contact your local government office.
All other accessories should also be disposed of in accordance with your local waste classification guidelines.

About the Battery
● If the operating time becomes extremely short even after charging, the battery may be at the end of its life. Do not use batteries with deteriorated performance. (Battery life may vary depending on usage and storage conditions. The estimated number of charge/discharge cycles is 300 to 500.)
● The battery is not replaceable.

For safety precautions regarding lithium-ion battery products, please refer to the following:
- Do not immerse in water or use as a power source.
- When disposing of this product, please follow your local municipality's waste separation rules.

The product (main unit) contains a built-in lithium-ion battery. When disposing of this product, please follow your local municipality's waste separation rules.

How to Store

● When not in use, turn off the power and untwist the USB cable.
● Do not store the product near fire, under direct sunlight, in bathrooms, inside cars, or in other high-temperature or high-humidity environments. Store it in a clean, well-ventilated, shaded area.
● Keep the product out of reach of children.
● If you do not plan to use the product for an extended period, fully charge the main unit before storage. Operate and recharge the product about once a month. If left uncharged for a long time, the battery may become unable to charge. Due to the nature of lithium-ion batteries, a gradual decrease in charge will occur even during storage.

● To prevent dust accumulation during long-term storage, cover the product with a cloth or similar material.
● You may roll it into a compact shape and fasten it with the hook-and-loop fastener for storage. (Figure 1)



Figure 1

CAUTION

Do not wrap the USB cable around the main unit during storage.
This may put stress on the cable, causing wire breakage, short circuits, or damage to the main unit or accessories.

Do not use any AC adapter other than one with a DC5V/1A output and Class II electric shock protection (□).
This may cause malfunction or accidents.

At the beginning of use, start with weak intensity of the pressure for about 5 minutes.
This may cause injury or health issues.

Limit each session to no more than 10 minutes, and the total daily usage to no more than 30 minutes.
Excessive use may be counterproductive or cause health issues. Avoid continuous use by multiple people.

CAUTION

If the heating function feels too hot, turn it off.
This may cause burns. Individuals who have reduced temperature sensation or are sensitive to heat should exercise extra caution.

When using a power tool, be sure not to exceed the rated capacity.
This may cause fire or damage.

Unplug the USB cable from the outlet when not charging.
Deterioration of electrical parts may cause electric shock, leakage, or fire.

If a malfunction or issue occurs, contact the retailer or our customer support.
Do not attempt disassembly or repair yourself. Such actions are not covered under warranty.

How to Clean

● Always clean the product only after it has completely stopped operating.
● Do not wash the main unit or accessories with water (excluding the detachable cover).
● To clean the main unit or accessories (excluding the detachable cover), wipe with a dry, soft cloth. If dirt is noticeable, soak a soft cloth in a solution of lukewarm water and mild neutral detergent, wring it out thoroughly, and gently wipe the surface of the main unit (fabric part).
● Do not use solvents such as alcohol, thinner, or benzene. This may cause malfunctions, cracking of parts, or discoloration.
● Only the detachable cover is washable. Follow the washing instructions below.

Care Instructions for the Detachable Cover
● Wash separately from other clothing and use a laundry net. ● Use a mild neutral detergent that does not contain bleach or fluorescent whitening agents. ● Do not leave it wet for an extended period, as this may cause color transfer or discoloration. ● After washing, hang it to dry in the shade while reshaping. ● Do not tumble dry. ● The cover may fade, so handle with care.

How to Use

3 Press the air level adjustment button () to adjust the air intensity.
The air level can be selected from the level.
*Each press of the button changes the level in the sequence of "Level 1 → Level 2 → Level 3 → Level 4 → Level 5 → Level 1 → Level 2."
*Start with a lower level and increase the air intensity gradually once you become accustomed to the product.

4 To use the heating function, press the heat level adjustment button ().
The heat level can be selected from 3 levels.
*Each press of the button changes the level in the sequence of "Low → Medium → High → Heat off."
*If the fit feels too tight, loosen it slightly and if the intensity of the pressure is too weak, select a bit more firm; this product can adjust the fit.

5 To turn off the power during use, press and hold the power button ().
The power will automatically turn off 10 minutes after the start of use.


Package Contents

If any items are missing from the package contents, please contact our customer support.

- Main Unit (MYTREX RAKUNO FIT) 1 unit
- Detachable Cover** 1 piece
- Charging Cable 1 piece
- Instruction Manual 1 copy

*1. The detachable cover is pre-attached at the time of shipment.
The detachable cover is washable. Please refer to "How to Clean" for details.
**The illustrations in this manual are for reference only and may differ from the actual product.
*Design and specifications are subject to change without notice.
*This product is not a medical device.

Names of Parts

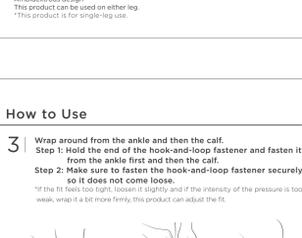

How to Use

1 Press and hold the power button () to turn on. The power will turn on and the indicator lights on the control panel will illuminate.
*Default setting is Course 1, Heat off.

2 Press the power button () briefly to select a course.
The 3 types of courses can be selected. Each press of the button changes the course in the sequence of "Course 1 → Course 2 → Course 3 → Course 1 → Course 2 → Course 3."
*The indicator lights are displayed when the power is turned on.

Course Name	Care Description
Course 1: Self-heating course	A sensation of being gripped and squeezed by the calf firmly.
Course 2: Pressure course	A sensation of being gripped by both hands and pressed for them.
Course 3: Refresh course	The care involves repeated tightening and releasing of the calf.

LED COLORS	Heat level indicator	Course indicator	Air level indicator
Off	Heat off	Course 1	Light blue Level 1
Yellow	Low	Course 2	Blue Level 2
Orange	Medium	Course 3	Green Level 3
Red	High		Yellow Level 4
			Red Level 5

Heat level adjustment Battery level indicators Air level adjustment button
Power/Course selection button
Press and hold: Turn the power on/off
Press briefly: Selects the course

Before Use (How to Charge)

When using the product for the first time or when the battery level is low (indicated by LED) following use of the battery level indicators, turn off the main unit before charging.

● Estimated Operating Time* ● Battery Level Indicators

Charging time	Operating time	LED 1	LED 2	LED 3	LED 4
Approx. 4 hours	Approx. 5.5 hours	●	●	●	●

*Operating time is based on usage under normal conditions (for Level 1) or low heat. When used under "Course Vibration (High heat)", the operating time is approximately 1.5 hours (at 25°C in a shaded environment). Charging and operating times may vary depending on battery level, charging conditions, and usage environment. Operating time may shorten in colder environments (below 10°C).
**Recommended charging temperature is 15°C to 35°C. Charging outside this range may result in longer charging times or charging failure.
**If charging for the first time, or after a long period of non-use, battery runtime may be shorter. Charging and usage multiple times will restore normal function.

How to Attach the Detachable Cover

The detachable cover can be removed from the main unit and is attachable to the main unit before use.

*For care instructions of the detachable cover, please refer to "How to Clean" (page 26).
*The detachable cover is pre-attached at the time of purchase.



How to Use

CAUTION
*Before use, check that the detachable cover fabric is not torn. Do not use the product if there is any damage.
*If using for the first time or after a long period of non-use, check that the buttons on the control panel and the timer (10 minutes) are functioning properly.
*Masks may remain on the skin after use. If this concerns you, wear socks or pants and then use the product over them.
*The barefoot can be used on either leg.
*This product is for single-leg use.

1 Connect the Type-C connector side of the included USB charging cable to the charging port of the main unit.

2 Connect the Type-A connector side of the included USB charging cable to the AC adapter (sold separately).
*Use a commercially available AC adapter that is able to output 5V/1A electric shock protection (□). Also, use an adapter with a rated output of DC5V/1A.
*Charging via a computer's USB port is possible, but some computers or mobile batteries may not support charging. If charging fails, try a different power source. Do not use a USB fast charger as it may cause malfunction.



Before Use (How to Charge)

3 Plug the AC adapter into a wall outlet.
● Charging begins, and LED1 to LED4 will blink in sequence, then light up. If any of the LEDs do not light up, charging is not being performed correctly. Please check if the connections are correct according to the procedures in "Before Use (How to Charge)."
*The product cannot be used while charging.

4 When charging is complete (all LEDs to LED4 light up), unplug the AC adapter from the outlet.

How to Attach the Detachable Cover

The detachable cover can be removed from the main unit and is attachable to the main unit before use.

*For care instructions of the detachable cover, please refer to "How to Clean" (page 26).
*The detachable cover is pre-attached at the time of purchase.

